
PETER



**BAUSCH**

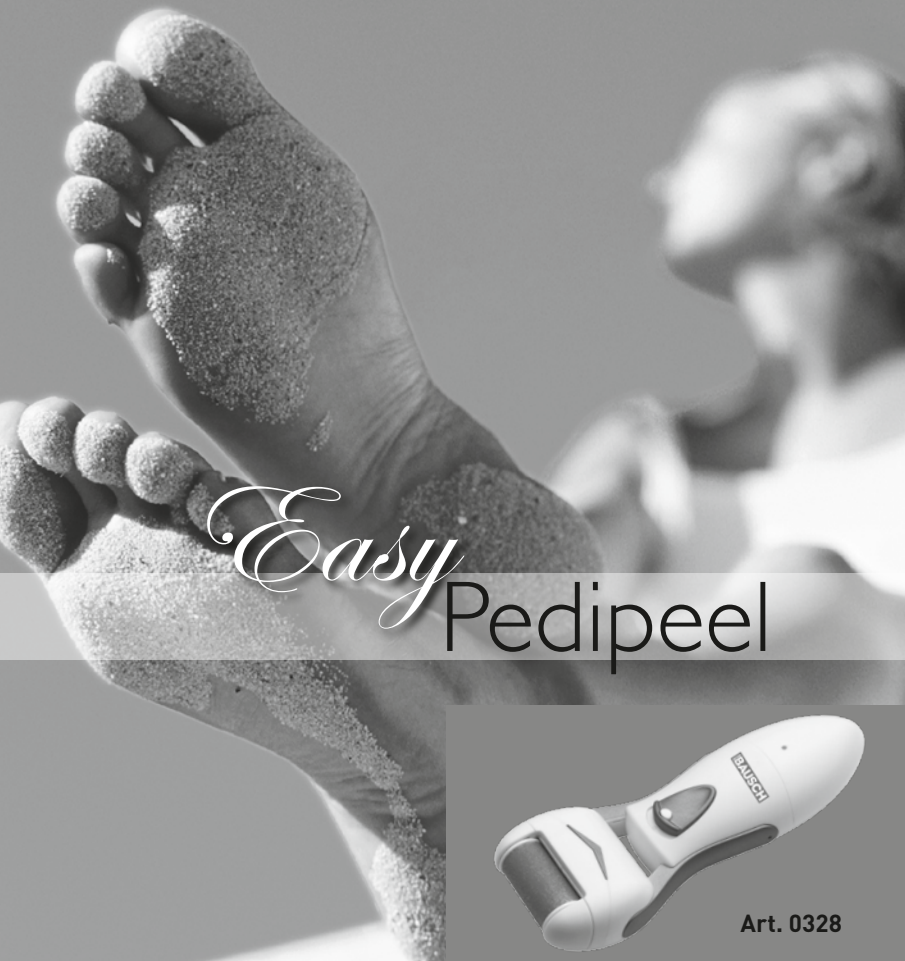*Easy* Pedipeel

Art. 0328

Gebrauchsanweisung • Instructions for use • Mode d'emploi  
Istruzioni per l'uso • Instrucciones para el uso • Gebruiksaanwijzing

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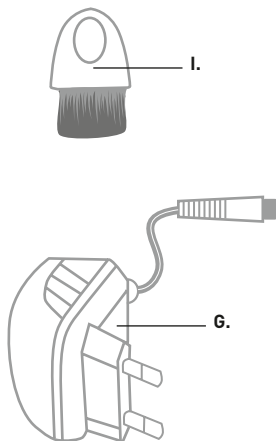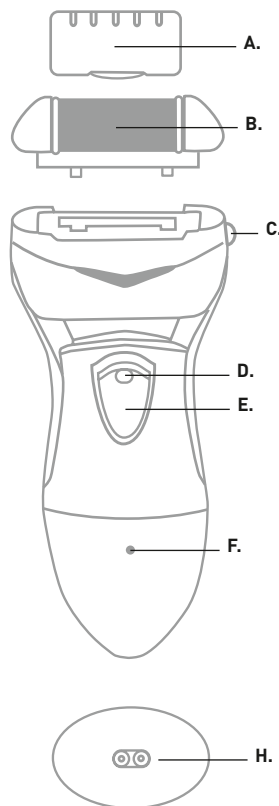- A. Protective cover
- B. Turning grinding head attachment (rough + fine)
- C. Release button, for changing the grinding head attachment
- D. Safety lock
- E. On /off switch
- F. Charging indicator lamp
- G. Mains adapter
- H. Connection point for charging cord
- I. Cleaning brush

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gb

Dear customer,  
thank you for purchasing our electronically powered device for a perfect pedicure! This manageable callus grinder EasyPedi Peel is especially made for removing dry, rough and also thickened callus softly. Due to its turning grinding head (40 revolutions a second) you are able to get a gentle, yet effective pedicure at an angle of 300°. This device is the perfect alternative to conventional methods about removing callus. It's chargeable, easy to clean, well-fitting for every hand and usable as easy as pie!

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Please read all available documentation (instruction of use and safety instructions) carefully before using this product! We assume neither responsibility nor liability for accidents, damages and malfunctions resulting from inobservance of these guidelines. The instructions have to be stored in a safe place and have to accompany the product at all times, especially if product gets passed on to a third party!

## SAFETY INSTRUCTIONS

Danger! The device is not to be used while taking a bath or a shower. Also, do not use the device over a fully filled sink! There's a risk of electrical shock!

- To be used under the direct supervision of an adult! This device is not a toy! Injury and electrical shock hazard!
- Devices with the nominal voltage 220/230V, 50Hz are only allowed to be plugged in appropriate isolated ground receptacles. Device may overheat!
- For cleaning, never immerse the device in water. There's a risk of electrical shock due to penetrated water!
- For operating only use the included mains adapter. Fire hazard!
- Do not bend or crush the connecting cable nor connect it to heated areas. Danger of electrical shock due to damaged cable insulation!
- Check the device every time you use it – do not use the device if it appears to be damaged or faulty. Injury and electrical shock hazard!
- **Attention!** Do not open nor damage by forcefully impacts the device. Injury and electrical shock hazard!

## INITIAL OPERATION

The EasyPedi-peel was developed to remove rough and thickened skin on feet and balls of the foot softly. Do not use it on other parts of the body!

1. Before first time usage or after using the device after a long period of non-operation make sure to fully charge the device (app. charging time 10h). Connect the tail end of the the charging cable (G) to the connection point (H) at the device. Due to the plugs special shape there is only one way to insert it correctly (+ in + and - in -)! Now plug the mains adapter into the receptacle, the charging indicator lamp provides information about the charging status.
  - You are able to use the EasyPedi-peel for 30 minutes if the battery is fully charged.
  - Regularly charge the device for 10h but never more than 16h.

2. Wash and clean your skin and dry it.

**Attention:** Never use the device on torn, ripped or bleeding skin.

3. Be sure that the turning grinding head attachment (B) is inserted correctly and engaged.
4. Remove the protective cover (A) from the attachment (B).
5. Press the safety lock (D) and simultaneously push the on/off switch (E) forward to turn the power on.
6. Circle the EasyPedi Peel over the affected areas for 2-3 seconds. Do not press too firmly; otherwise you might dislike the feeling and the device will stop! Repeat this step until you have reached the softness you want to achieve.

**Attention:** Never circle over the same area longer than 2-3 seconds! Stop the process if your skin gets too thin or hot!

7. Turn off the device after finishing your pedicure.
8. Clean your skin from all remaining skin particles and use enriched foot lotion for soft and healthy feet.

#### ■ EXCHANGE OF THE ATTACHMENTS

The turning grinding heads are interchangeable and can be exchanged on demand.

1. The device has to be turned off.
2. Push the release button (C) on the side of the EasyPedi Peel and remove the grinding head (B) easily.
3. Attach another grinding head to the device. Make sure that the attachment is fully engaged. **This is just possible if installed properly!**

**TIP:** The small gear-wheel on the bottom side of the attachment has to point into the same direction as the gear-wheel of the top side of the device! If installed correctly, you will hear a "click" sound.

4. Place the protective cover (A) over the attachment.

**Tip:** Do not turn on the device if there is no attachment attached to it!

#### ■ CLEANING, CARE AND MAINTENANCE

**Tip:** To ensure a constant performance of the device, clean it after every usage thoroughly. You can rinse the attachment under running water without any concerns. Do not immerse the device into water! Do not clean it with chemicals! The device is maintenance-free, be careful with it!

1. The device has to be turned-off during cleaning.
2. Push the release button (C) on the side of the EasyPedi Peel and remove the grinding head (B) easily.
3. Remove all remaining skin particles with the cleaning brush (I). You are able to clean dust and skin particles from handle and attachment with it.
4. Rinse the head piece of the device with water. While doing so turn the roll. You can also use soap for cleaning the grinding head attachment. You may disinfect the attachment by using disinfectant spray on it.
5. Dry all the parts with a towel.
6. Attach the grinding head to the device. Make sure that the attachment is fully engaged. **This is just possible if installed properly!**
7. Place the protective cover (A) over the attachment.

**Tip:** The small gear-wheel on the bottom side of the attachment has to point into the same direction as the gear-wheel of the top side of the device! If installed correctly, you will hear a "click" sound.

**Tip:** Do not turn on the device if there is no attachment attached to it.

#### ■ REPAIR INSTRUCTIONS

**Attention!** Electric devices are only fixed by qualified personnel. Improperly repairing can result in severe consequential damage.

#### ■ DISPOSAL

This product as well as any other electronic and electrical appliances, which have reached the end of their lifetime, must be collected and disposed of separately from household waste material.

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