
PETER **BAUSCH**

*Professional*  
Manicure & Pedicure

Art. 0379P

Gebrauchsanweisung • Instructions for use • Mode d'emploi  
Istruzioni per l'uso • Instrucciones para el uso • Gebruiksaanwijzing

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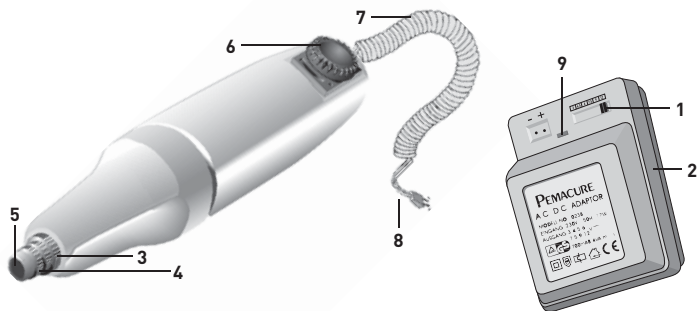1. Speed regulator
2. Adapter
3. Knurled screw
4. Nut
5. Location hole
6. On/Off (Speed regulator)
7. Cable
8. Small plug
9. Pilot lamp

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Dear customer,  
congratulations for having purchased an electric manicure and pedicure set for perfect hands and feet and thank you for having chosen our product, surely due to the quality and the design as these are our principal aims together with the satisfaction and enjoyment of our clients.

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Please, read and follow these instructions of safety before using this appliance.

Our warranty doesn't cover damages and/or injuries caused by non-observance of these safety norms and instructions. Keep this booklet of instructions in a safe place and adjoint it to the appliance if or when you should give it away.

#### STARTING UP

Hold the knurled screw [3] with one hand. Unscrew the nut [4], with your other hand, turning it in a clockwise direction. Insert the attachment into the location hole [5]. To tighten the nut [4] turn in an anticlockwise direction.

Connect the appliance to the adapter [2] with the small plug [8], conforming the poles (+/-). Plug in the adapter. The pilot lamp [9] indicates proper connexion to mains. Adjust the speed of rotation on the adapter for the best performance of the appliance on the adapter according to the result wanted (highest voltage on 9 or 12 Volt).

**Select the speed of rotation (on the adapter), low for sensitive skin:** Turn the On/Off regulating wheel to the right past the click-in point. The machine will start up at a low speed. Turn the On/Off regulating wheel further to the right in order to steplessly raise the speed. The highest speed is on the right-hand stop of the regulating wheel.

#### USE OF THE ACCESSORIES

Hands and feet can be treated a lot easier after a warm bath (with or without additives). Nails are then more flexible, and hard skin is easier to remove.

**Pointed abrasive tip:** for shortening and shaping finger and toe-nails.

**Cylindrical abrasive tip:** to remove hard skin from fingertips, toes, heels and sole.

**Cleaning and polishing brush:** use this to remove skin and nail dust, or to polish the nail top. Will also clean underneath nail-ends.

**Felt polishing tip:** to polish the nails prior to application of lacquer. Using regularly this polisher will prevent nails to split or break.

**Diamond cutter:** to remove hard skin from fingertips, toes, heels and sole.

**Diamond cutter flame shaped:** for the removal of horny skin or corns and partridge eyes.

**Diamond cutter:** for dealing with ingrained nails as well as for removing the not very aesthetic cuticles around the nails.

#### IMPORTANT SAFEGUARDS

Never use the appliance in the bath, or near a water-filled hand-basin.

- Do not let children use the equipment without supervision.
- Connect the adapter exclusively to AC mains of 230 Volts ~ 50Hz.
- Never immerse the manicure/pedicure appliance in water.
- Connect the appliance exclusively to the adapter supplied with it.
- Never bend or squeeze the cable, nor let it get in contact with hot things.
- After each use unplug the adapter from the mains.
- Before each use check the appliance, in case of damage of any kind avoid using it.
- **IMPORTANT:** The appliance cannot be opened, not even by twisting it.

#### CLEANING AND MAINTENANCE

For cleaning just wipe the manicure/pedicure appliance with a dry, soft cloth. Don't ever put the set in water or use detergents. Except for the polishing felt and the brush, you can use water, soap and a brush to clean the other accessories. After the cleaning it is advisable to spray on some disinfectant. Let the tips dry in order to avoid rusting. The set needs no maintenance, just take good care of it.

#### REPAIRS

Caution! Electrical appliances should only be repaired by qualified technicians, as otherwise serious damage may be caused.

#### WASTE DISPOSAL

At the end of its working life the appliance must not be disposed of as urban waste. It must be taken to a special local authority differential waste collection centre for electrical and electronic products.

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